




P05 Portable Oxygen concentrator

Make life further and more possible




Portable Oxygen Concentrator P05

The P05 portable oxygen concentrator is a device to provide oxygen therapy, using pressure swing adsorption (PSA) technology and highly efficient molecular sieve, delivers up to 96% pure oxygen to people who require supplemental oxygen. The P05 provides oxygen in pulse mode, which improves oxygen utilization efficiency.

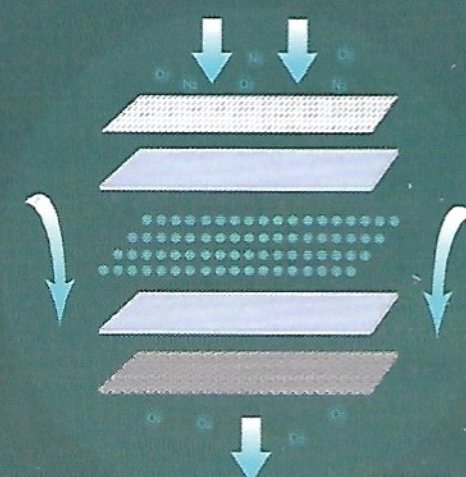


Oxygen enrichment output

- ▶ 1-5 settings, can be adjusted to meet different needs
- ▶ Replaceable dual-cylinder molecular sieve outputs high-purity medical-grade oxygen
- ▶ Stable output of $93 \pm 3\%$ oxygen-enriched air
- ▶ Oxygen is provided when the user inhales, reducing oxygen waste
- ▶ With a sensitivity of $0.2 \text{ cmH}_2\text{O}$, the pulse oxygen supply is more synchronized with breathing

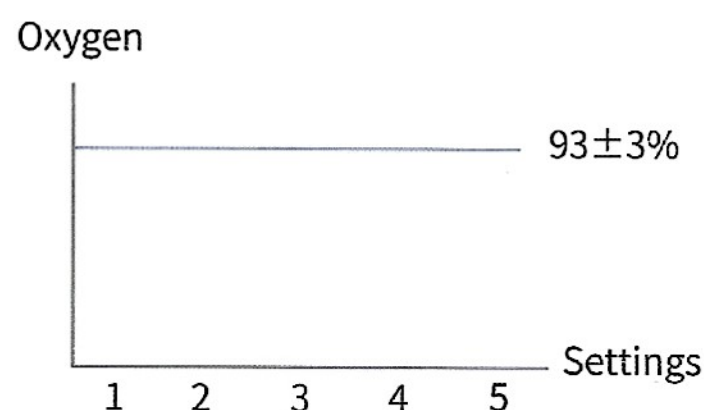

Backup continuous mode

- ▶ Providing backup pulse mode and **backup continuous mode**
- ▶ In the backup continuous mode, if no spontaneous breathing is detected for 30 seconds, the backup continuous mode will be triggered, supplying oxygen at a flow rate of 1000ml/min
- ▶ During the backup continuous oxygen supply, if 3 breaths are detected within 2 mins, the device would switch to pulse mode oxygen delivery

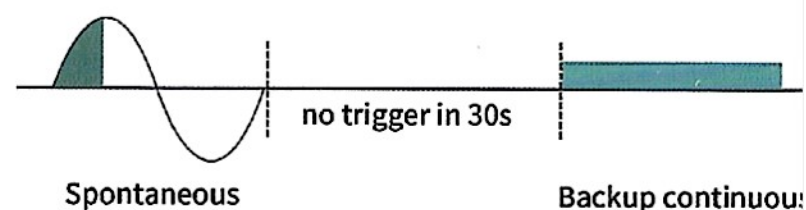

Smart for outdoor

■ Light weight design

Less than 2.0kg

Only 20.8cm tall

P05 is not only suitable at home and indoor environment, but also for vehicle or other outdoor mobile environment.

■ Long-lasting duration

The detachable lithium battery can be charged independently through a type-C adapter.

Two capacities of batteries: BAT-01 and BAT-02

BAT-01 has a battery life of 6 hours, and BAT-02 can last up to 12 hours, providing greater peace of mind when going out.

■ Easy operation

2.8" color screen with one-button power on/off switch



Caring guardianship

- ▶ 12 graded alarms with real-time monitoring and indication through synchronized display of indicator lights
- ▶ The alarm log function records alarm events for later review

Alarm log

No breath	22-12-12 13:30
Low battery	22-12-12 13:30
Regular Maintenance	22-12-12 13:30
Internal failure	22-12-12 13:30
High Pressure	22-12-12 13:30

Specifications Sheet

Model	PO5			
Pulse volume	1	220ml (±15%)	4	880ml(±15%)
	2	440ml(±15%)	5	1000ml(±15%)
	3	660ml(±15%)		
Dimension	188mmx98mmx208mm			
Weight	2.0kg (with battery)			
Screen	2.8"			
Mode	Pulse, backup continuous mode			
Oxygen concentration	93±3%			
Warm-up	2 minutes			
Battery duration	BAT-01: 1	360mins	BAT-02: 1	720mins
	2	195mins	2	390mins
	3	170mins	3	340mins
	4	150mins	4	300mins
	5	120mins	5	240mins
Charge time	BAT-01:>2.5 hours		BAT-02:>5hours	
Operation conditions	Temperature 5°C-40°C Relative humidity: 35%-93% (non-condensing)			
Storage conditions	Temperature -20°C-55°C Relative humidity: 10%-93% (non-condensing)			

SHENYANG RMS MEDICAL TECH CO., LTD

ADD:NO. 10-4, JINHUI STREET, HUNNAN DISTRICT, SHENYANG

TEL:024-31682686

E-MAIL:cooperation@rms-care.com